

Certificate of Analysis

Eurolab Peptides

Eurolab Peptides
Belgium

Liquilabs s.r.o.

Inovační 122
252 41 Zlatníky-Hodkovice
Czechia
www.liquilabs.cz



[Verify Results Online](#)

Sample Identification

Sample Name Epithalon 10 mg
Batch Number
Date Published 2026-07-17 17:35

Results for LYO-0350

Peptides	Result	Unit	Uncertainty	Acceptable Range
Epithalon Assay Peptide Screening 0.1% TFA	6.03	mg	[± 0.03]	
Epithalon Purity Peptide Screening 0.1% TFA	99.2	%	[± 0.5]	
Epithalon Identification by Spectrum (FTIR) Peptide Screening 0.1% TFA	981		[± 5]	
Epithalon Identification by RT Peptide Screening 0.1% TFA	1.000		[± 0.005]	

	Method Specification		
Determination of identity, content and purity of Epithalon			
<i>Document number</i> EPI_006_2026	<i>Superseded document</i> -	<i>Number of pages</i> 4	

1. Content Assessment

1.1. Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

1.2. Chromatographic conditions

Chromatographic conditions	
Eluent A	0.1% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.1% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.4 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	60°C
Column	Phenomenex Biozen Peptide Polar C18, 150x2.1mm 3µm
Detection wavelength	225nm


Gradient Program		
Time [min]	A [%]	B [%]
3	100	0
18	45	55
20	1	99
21	1	99
21.01	100	0
26	end	

1.3. Sample preparation

Whole amount of container was dissolved in 2mL of water (LCMS Grade). 100 µL of sample was transferred to HPLC vial and diluted by 900 µL water (LCMS Grade) and submitted for analysis.

1.4. Calibration curve

Calibration curve detail	
Quantitative method	External Standard
Calibration Type	Linear
Number of calibration points	5
Force through Zero	Enabled
Weighting Method	None

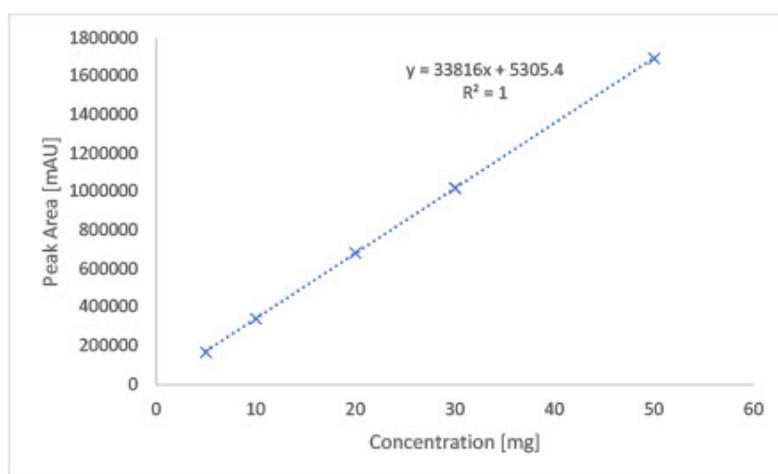

2. Purity assessment

2.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

2.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.1% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.1% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.4 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	60°C
Column	Phenomenex Biozen Peptide Polar C18, 150x2.1mm 3µm
Detection wavelength	225nm

Gradient Program		
Time [min]	A [%]	B [%]
3	100	0
18	45	55
20	1	99
21	1	99
21.01	100	0
26	end	

2.3 Purity assesment

Purity of compound assesed by area normalization method, comparing area of each peak to sum of area of all peaks detected at wavelenght of 225 nm.

Analysis Report



LIQUILABS

Analytical services


Analysis of quantity and purity of active ingredient by UHPLC with UV detection

Sample Information


Injection Volume : 2
Data File : LYO-0350__003.lcd
Method File : Peptide screening polar Group C.lcm
Date Acquired : 7/15/2026 12:29:54 PM

Chromatogram

uAU

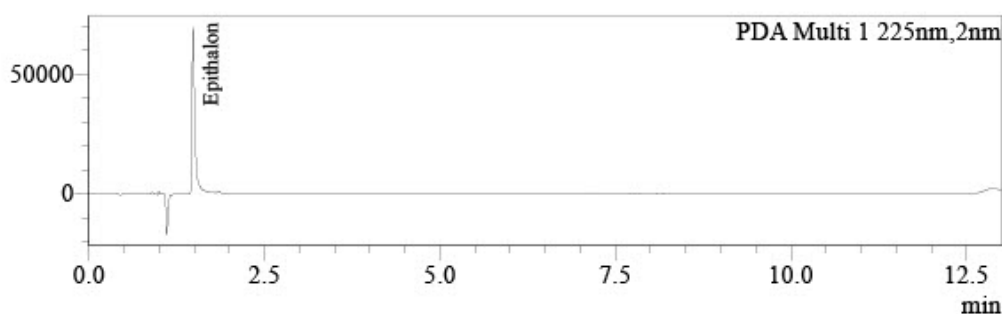

uAU

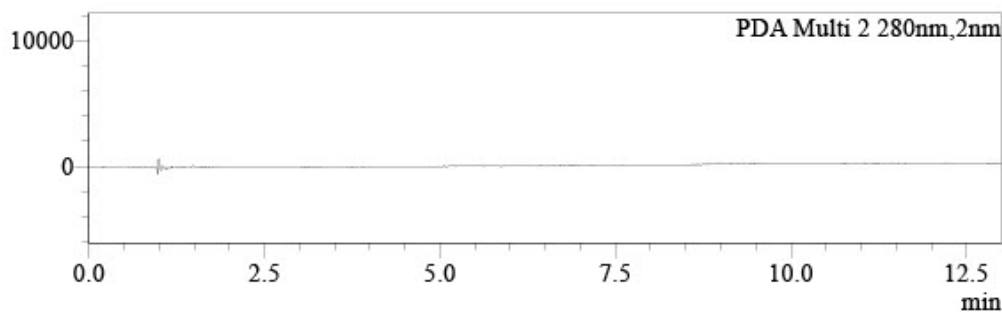

Peak Table

PDA Ch1 225nm

Name	Ret. Time	Area	Conc.	Unit	Area%
Epithalon	1.480	209052	6.025	mg	99.236
	1.843	1609	0.000		0.764
		210660			100.000

Peak Table

PDA Ch2 280nm

Name	Ret. Time	Area	Conc.	Unit

Responsibles



Mr. Ján Galbavý
CEO

Analysis results relate only to the samples tested.

This document shall not be reproduced except in full, without the written approval of Liquilabs s.r.o.